

6.4 kb

Page N _____

urpose To amplify 6.4 kb and 8.0 kb from plasmid
used F + R (non du) primers

510 µl Rx 200 µl dH₂O each cycling 94°; 1'
14 µl primers
2 mM Mg used buffer B (94°; 30" 60°; 45" 72°; 3") 25
Template ?
1 µl enzyme per mixed (1:0.01)

used enough premix for 20 Rx:

6.4 kb

all done in duplicate.

included purified }
prep at a known }
concentration }
con tried 50 µg + 100 µg
(Tag 50) just one.

miniprep, unknown }
concentration (from }
the amount of colonies }
in 1/100 dilution) }
used .5 µl and 1 µl
Con should be quite }
high in the miniprep }
diluted to 60 µl }

plasmid - picked a single isolated colony directly
into the reaction mix containing all the rest of the stuff
done in duplicate


8.0

no purified stuff available

min prep - }
unknown con }
lot of colonies }
from 1/100 → 25 µl }
dilution }
- 5 and 1 µl
(out of 60 µl from
1.5 ml culture)
plasmid - 2, one in each - done in duplicate

To Page No. _____

passed & Understood by me,

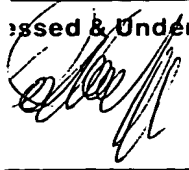

Date
1/9/95

Invented by


Recorded by

S. Sitaromun

Date

1/5/95

From Page No. _____

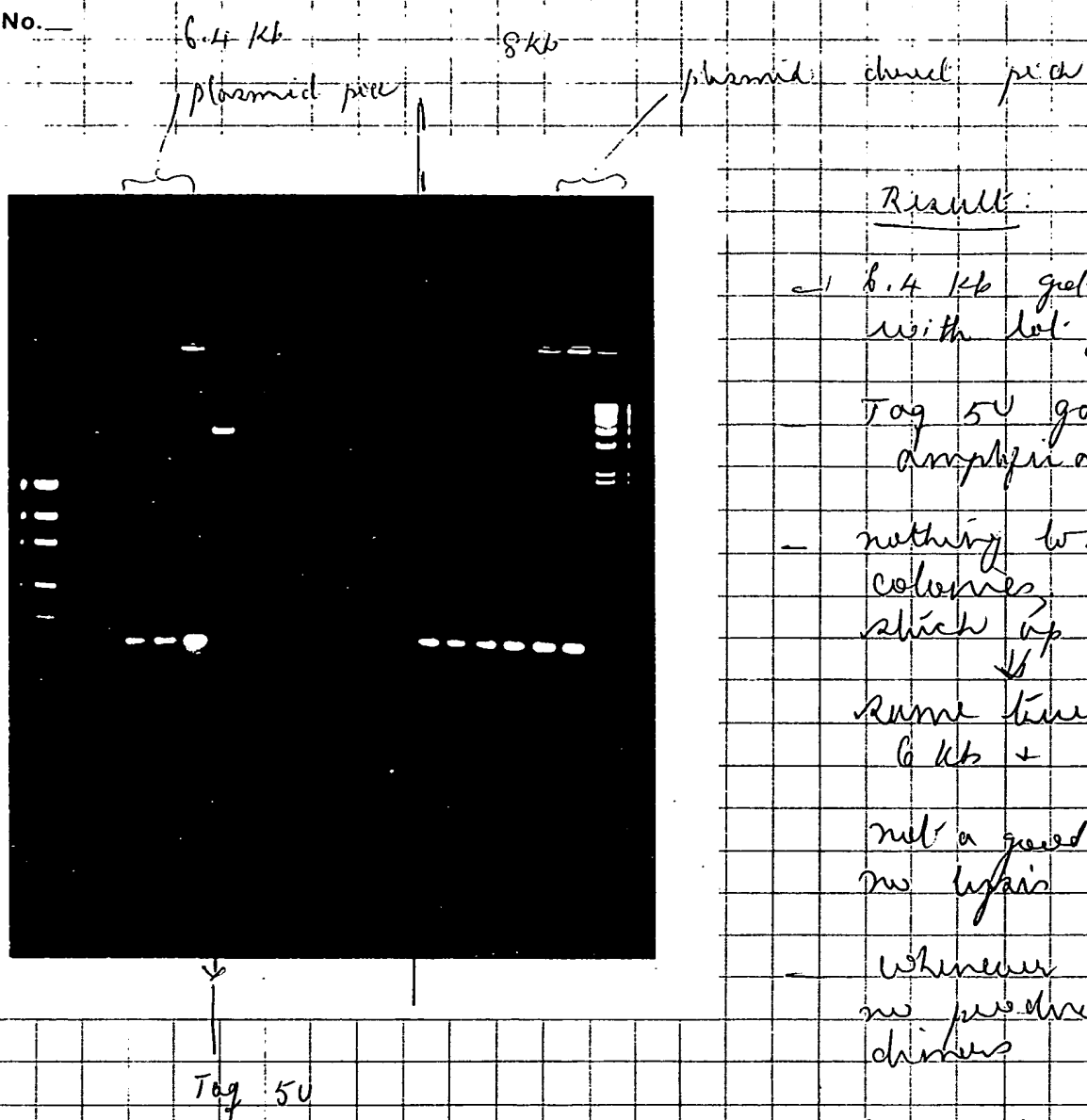Results:

- 6.4 Kb gel amplified with lot of mispriming

Tag 50 gave good amplification

- nothing to be seen for colonies, lot of stuff stuck up in the well

Same time for the 6.4 Kb + 8. Kb

not a good way to go, no lysis at all,

- whenever there was no product, lot of primers

- amount of primers to be found anyway.

* check alternate cycling conditions to get rid of mis priming

* lysis in PK and first reaction has been checked next

* make 6.4 Kb to work first

To Page No. _____

Witnessed & Understood by me, _____

Date _____

Invented by _____

Date _____

Recorded by _____

Dr. Subramaniam

1/9/25