

6.4 kb

Book N \_\_\_\_\_

153

Page N \_\_\_\_\_

urpose To amplify 6.4 kb and 8.0 kb from plasmid

used F + R (non du) primers

510 µl Rx 200 µl dH<sub>2</sub>O each

1.4 µl primers

2 mM Mg

Template ?

1 µl enzyme per mixed

used buffer B

cycling: 94°; 1'

$\left( \begin{array}{cc} 94^{\circ} & 30'' \\ 60^{\circ} & 45'' \\ 72^{\circ} & 3' \end{array} \right) \times 25$

used enough premix for 20 Rx:

6.4 kb

all done in duplicate.

included purified prep at a known concentration

con tried 50 pg + 100 pg

(Tag 50) just one.

miniprep, unknown concentration (from the amount colonies in 1/100 dilution)

Con should be quite high in the miniprep diluted to 60 µl

used .5 µl and 1 µl

plasmid - picked a single isolated colony directly into the reaction mix containing all the rest of the stuff done in duplicate

8.0

no purified stuff available

min prep

unknown con

lot of colonies

from 1/100 → 25 µl dilution

.5 and 1 µl

(out of 60 µl from 1.5 ml culture)

plasmid

2

one in each

done in duplicate

To Page No. \_\_\_\_\_

passed & Understood by me,

Date

1/9/95

Invented by

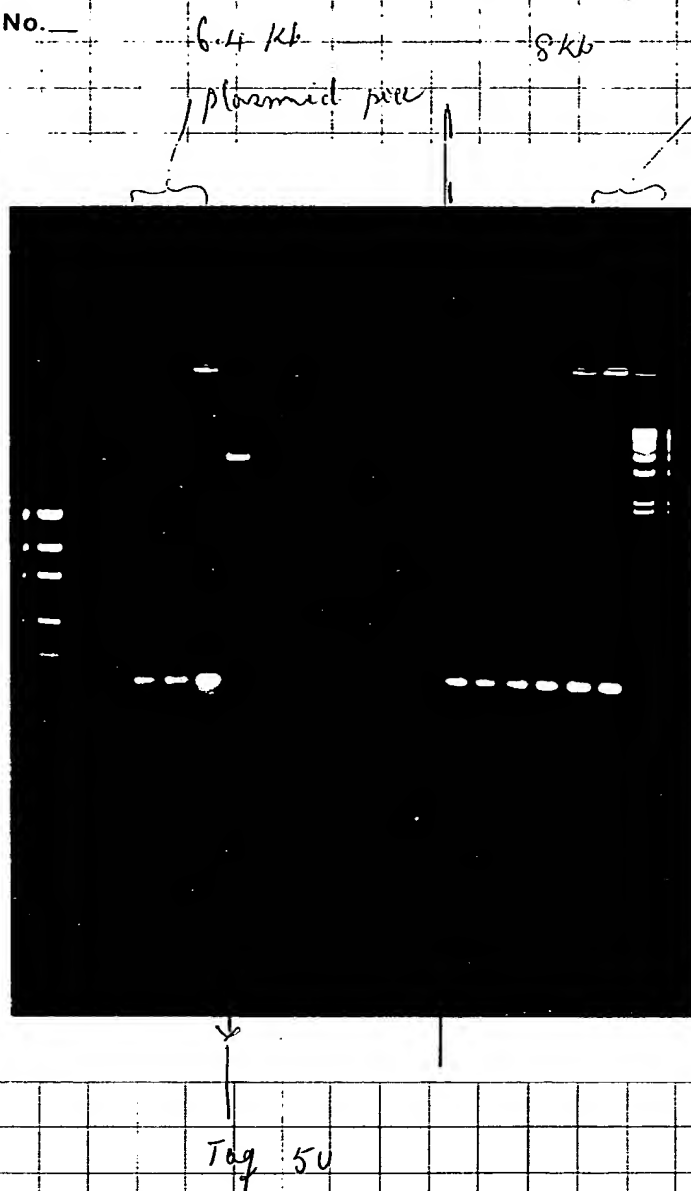
Recorded by

S. Sitaromun

Date

1/5/95

From Page No. \_\_\_\_\_

Result:

- 6.4 Kb gel amplified with lot of mispriming

Tag 50 gave good amplification

- nothing to be seen for colonies, lot of stuff stuck up in the well

Same time for the lot 6 Kb + 8 Kb

not a good way to go, no lysis at all,

- whenever there was no product, lot of primers

- amount of primers to be good enough.

\* check alternate cycling conditions to get rid of mispriming

\* lysis in PK and first reaction has been checked next

\* make 6.4 Kb to work first

To Page No.

Witnessed &amp; Understood by me,

Date

Invented by

Date

Recorded by

Dr. Subramaniam

1/9/25