

36

Project No. _____

Book No. _____

TITLE Clone both mutation in one.

From Page No. _____

First: Introduce both mutation in puc.

Digest (1) pucTne35 with BamHI → CIP → run gel

(2) pucTneFY with BamHI → Run gel




Purify largest frag. from pucTne35 and the smallest fragment from pucTneFY.

Ligate transform DH10B.

Plate 10% & 90% overnight at 30°C.

Grow: pTneFY 2L 1 Liter → 30°C 5 hrs.
Induce with IPTG (1mM) → Harvest.

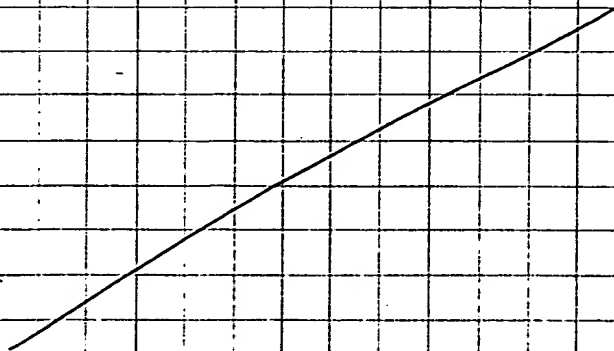

T Pag N

Witnessed & Understood by me,

[Signature]

Date

6/12/95

Invented by

Recorded by

[Signature]

Date

6/10/95