

T. neapolitana DLE (2779) SOM

Page N 48

October 13, 1994 (Thursday)

I infected DH12S cells with the ϕ from #5, #6, #7 (2ml cells grown in 2x YT + 5 μ l ϕ ; grown 37°C (air shaker) 16 hrs)

RF isolated by alkaline/SOS except 1 μ l RNase A (1mg/ml) added to prep at NH₄OH addition


RNA dissolved in 50 μ l T₁₀E₁

DIGEST SCHEME

(React 3)	HOP	10 μ l	✓	Incubated 37°C (heat block) 1 hour.
	10x Bfr	2	✓	
	DNA	7	✓	
(4 μ l)	Eco47III	1	✓	
	Total	20 μ l		


3% Agarose Gel (1xTAE); 190V/16

1. Run #2779 Eco47III DyeS

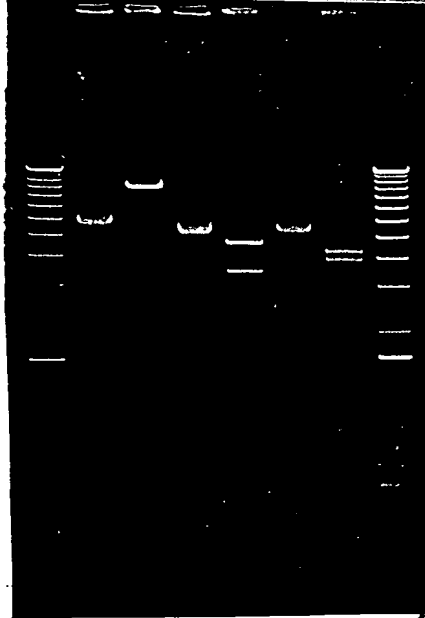

6 Kb
5 Kb
4 Kb
3 Kb
2 Kb

Comment:

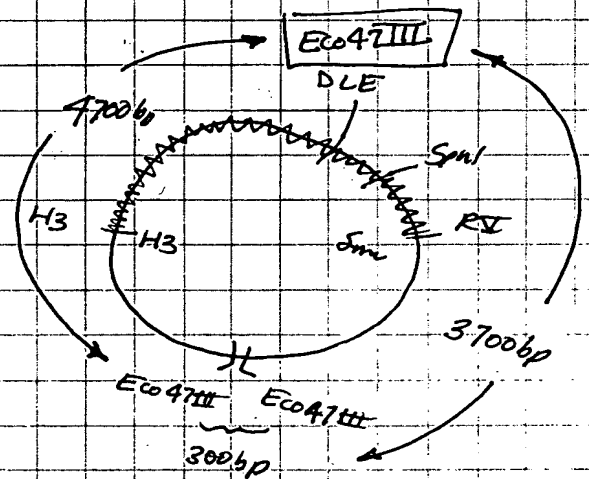

The bands are migrating at the expected distances for #6. There must have been an overabundance of some "component" causing the DNA to run faster. I will clone into the pUC vector.

To Page No. 50

ss d & Understo d by me,

Date

Invented by

Date

May forgo

10/24/94




Record d by

Pr j ct No. 20222

10-13-94

From Page No. 49

October 13, 1994 (Thurs)

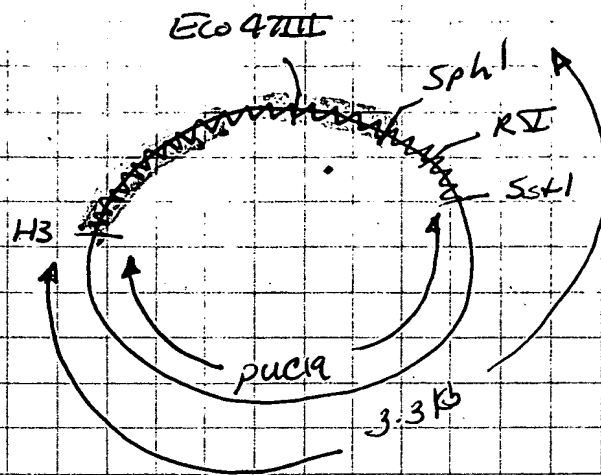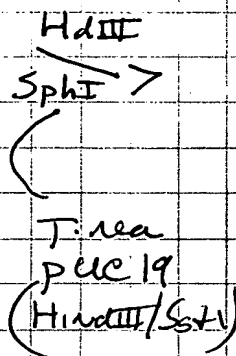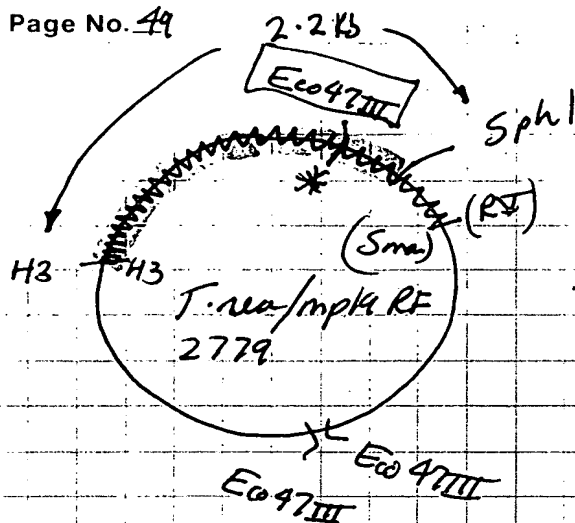


DIGEST SCHEME

		<i>T. nea/puc</i>	#6	<i>2779/mp19</i>		
(React 2)	HOM	12 μ e	✓	6 μ e	✓	Incubated 37°C (heat)
	10X Bfr	2	✓	2	✓	
	DNA	4	✓	12	✓	
(104 μ e)	HindIII	1	✓	1	✓	
(104 μ e)	SphI	1	✓	1	✓	
	Torac	20 μ e		20 μ e		1:10 \rightarrow 1:50

0.8% Agarose Gel (1XTAE), 190



10-24-94

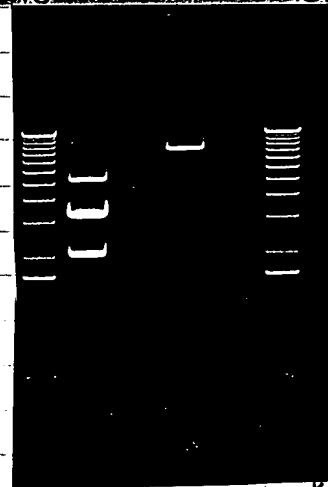

Comments: I should see a 3.3 Kb (desired fragment) and a 2.2 Kb fragment from the *T. nea/puc19* clone and I do. Unfortunately I should see a 7.9 Kb fragment and 2.2 Kb (desired fragment) from the *T. nea/mp19* (2779 #6) RF DNA and I don't. Both sites were present before I performed the mutagenesis (see p. 35) - I will have to repeat.

To Page No

Witnessed & Understood by me,

Date

Invented by

Date

M. Jones

10/24/94

Recorded by

D. J. Schmitt

10-13-94

ge N _____


January 25, 1995 (Wednesday)

I'm BACK!!

after a tour of duty with Joel's group, a trip to Aulba and a few week vacation off I am back and ready to administer the fatal blow to this project. I will finish sequencing this gene, mutagenize it to conform to our needs, and overexpress it so people can enough enzyme to swim in it and still have money left over for a cup of coffee and a copy of the New York Times!

How's that for an opening!

First things first. Let's reclone the region of the pol gene we are interested in mutagenizing. Rob and I have had no success with the last clone. Secondly, let's make the subclone more user friendly by eliminating the 3' BamHI site.

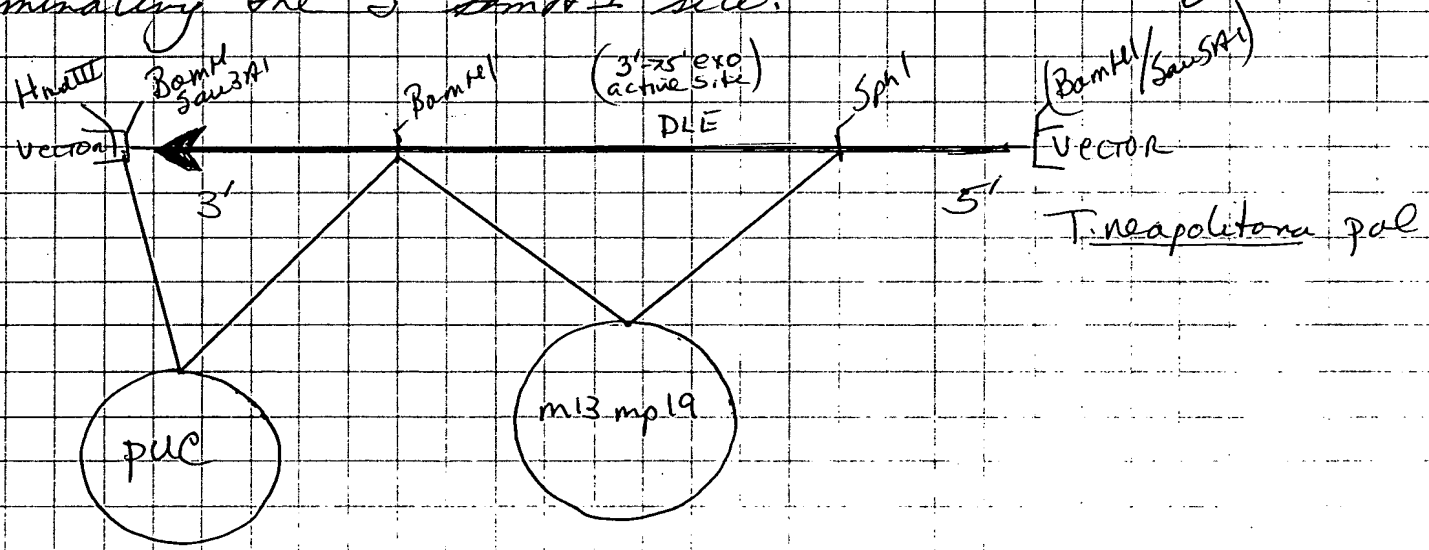

stoy BamHI/Sau3A
PCR
clone Hd3/BamHI

make ssDNA
D → A by SDM

To Page No. 52

ed & Understood by me, my forgo	Date 1/27/95 4/24/95	Inv nt d by Recorded by Schmidt	Date 1-25-95
------------------------------------	----------------------------	---------------------------------------	-----------------

From Page No. 51

January 25, 1995 (Wednes.)

T. neapolitana pSPORT DNA made by Michael Smith
(not the Nobel Laureate; the Hozobag)


DIGEST SCHEME

	<u>T. nea/pSPORT</u>	<u>mBmp19</u>	<u>pUC18</u>
(React 6) HOH	15 μ l ✓	13 μ l ✓	13 μ l ✓
10X B/R	2 ✓	2 ✓	2 ✓
DNA	1 ✓	3 ✓	3 ✓
(100 μ l) BamHI	1 ✓	1 ✓	1 ✓
(100 μ l) SphI	1 ✓	1 ✓	1 ✓
Form	20 μ l	20 μ l	20 μ l
(0.14 μ l) CAP			1 μ l

Incubated 37°C (heat-block) 1:00 \rightarrow 2:45

0.8% Agarose Gel (1XTAE)
190 Volts

BamHI/SphI
12/27/95

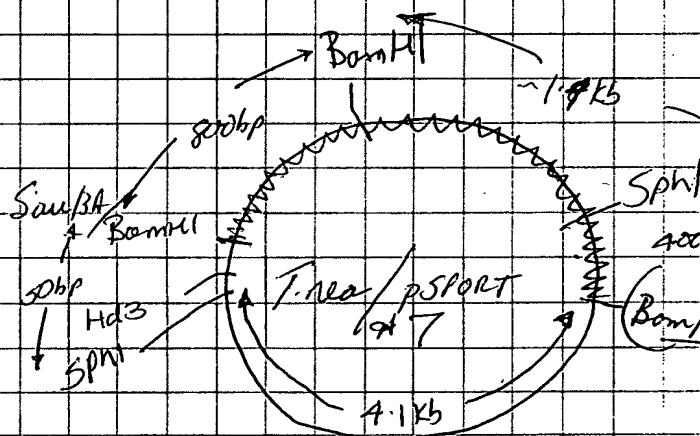

I forgot to run the 1Kb ladder.
what a hore bag!

Fragment T. nea/pSPORT 7 BamHI/SphI shows
the 4.5 Kb, 1.4 Kb, 0.8 Kb, 0.25 Kb.
Perhaps something partialled. Try again
but do it separately.

Witnessed & Understood by me,

Date

Invnt d by

Date

May Jones

1/27/95

R cord d by

Dr. Michael Smith

1-25-95

T. neapolitana SDM

ag N 52

January 26, 1995 (Thursday)

DIGEST SCHEME:

		<i>T. nea</i> /psPORT		m13 mp19 (~270ng/ul)	
1 (React 3)	HOH	24.5 μ l	✓	22.5 μ l	✓
	10X B/R	3	✓	3	✓
	DNA	1	✓	3	✓
F107 (104 μ l)	BamHI	1.5	✓	1.5	✓
	TOTAL	30 μ l		30 μ l	


Incubated 37°C (heat block) 8:04 → 9:08

3 μ l removed for analytical gel.


		<i>T. nea</i> /psPORT	mp19
DIGEST	27 μ l	✓	✓
1M KCl	2	✓	✓
HOH	9	✓	✓
1M SpH	2	✓	✓
TOTAL	40 μ l		

Incubated 37°C (heat block)
9:17 → 10:25

Agarose Gel (1XTBE)
190V 16

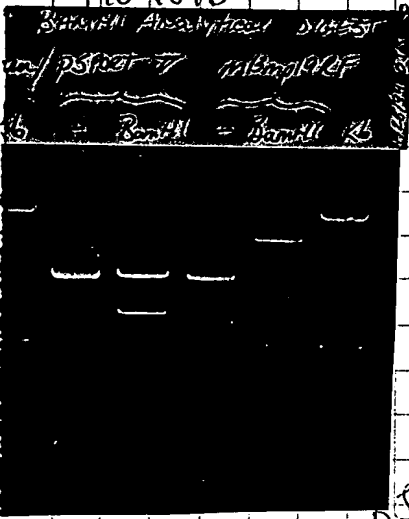

Comments

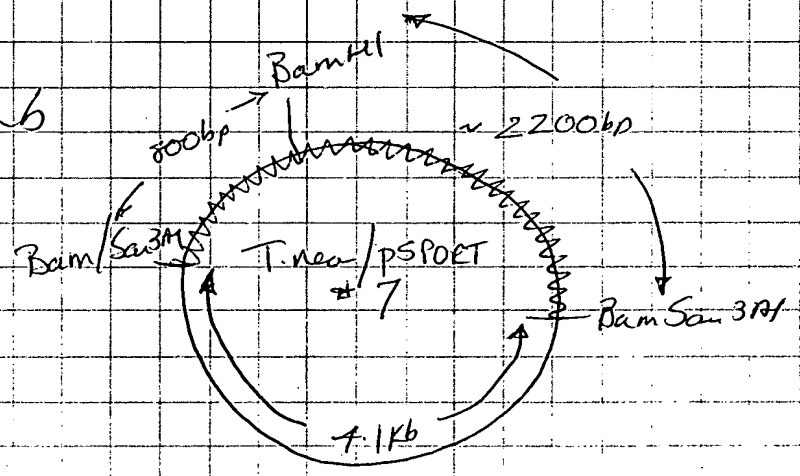

To Page No. _____

is d & Understood by me,

Date

Invent d by

Dat


Mer Torzo

1/27/95

Recorded by

Drumf. Schmidt

1-26-95

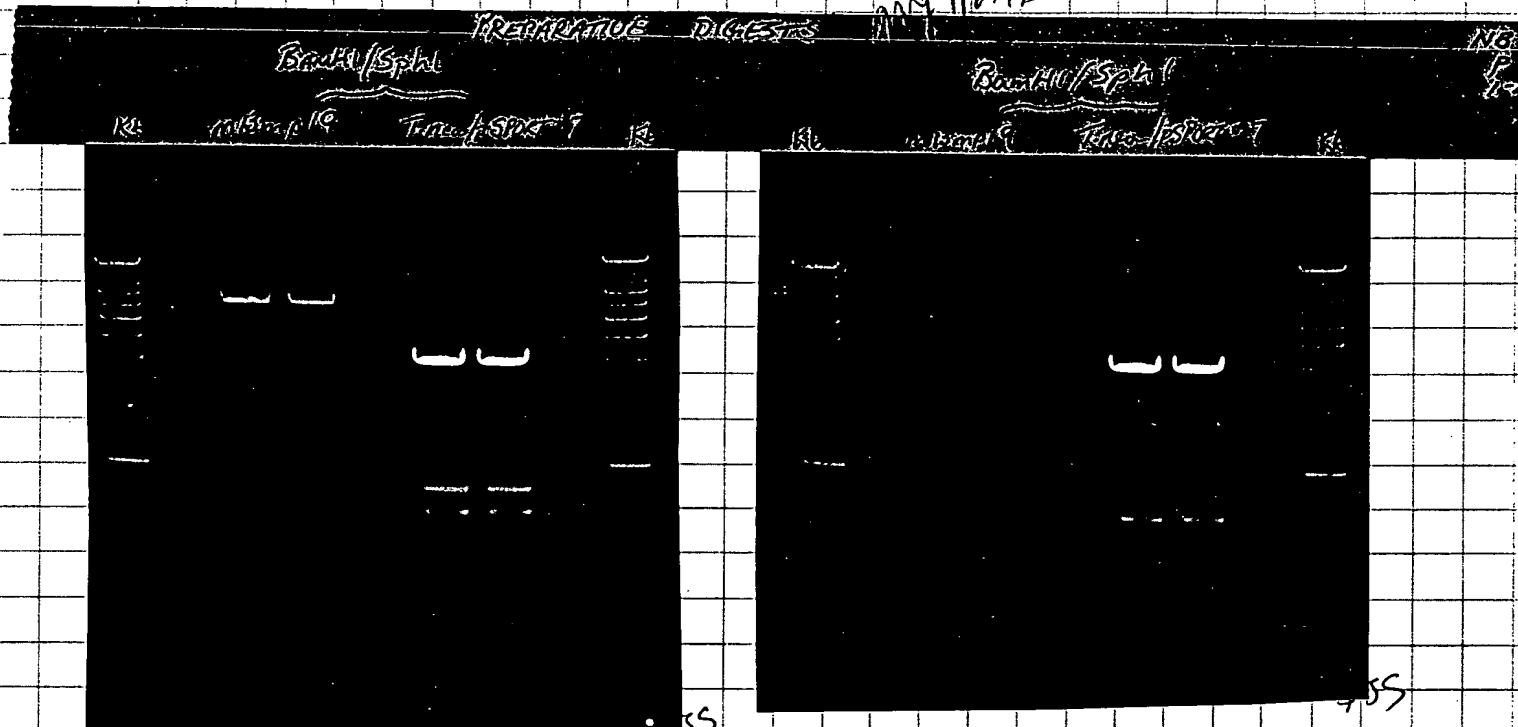From Page No. 53January 26, 1995 (th0.8% Agarose Gel (1X TAE); Run at 190 Volts1/27/95

Bands extracted from the gel and placed in the same tube. The DNA was purified away from the agarose using Gene Clean as described by the manufacturer (BIO-101)

DNA eluted in 14 μ l H₂O

LIGATION SCHEME

ETG 402 (Ligase)	DNA	14 μ l	✓
	5X Bfr	4 μ l	✓
	10 μ l (Ligase)	2 μ l	✓
	Form	20 μ l	

Incubated 22°C (room-temp)
2:15 → 3:15

→ 1 μ l ligation / 2 μ l for transformation

To Page

Witnessed & Understood by me,

Mary Loup

Date

1/27/95

Investigated by

Recorded by

David J. Kimmel

Date

1-26-95

T. neapolitana 50M


Pr j ct N .
B k N .

age No. 54

January 26, 1995 (Thursday)

DH10B Electrocompetent
20 µl DH10B Electrocompetent Cells + 1 µl (of a 1/3 dilution; See p. 54)
2.5 KV
1 ml LB, 37°C au shaker 20 min
↳ 10% applied to LB plate in 4 ml Soft Agar (0.7%)
90% + IPTG (1 mM) and X-gal 100 µl of 4%

incubated 37°C incubator

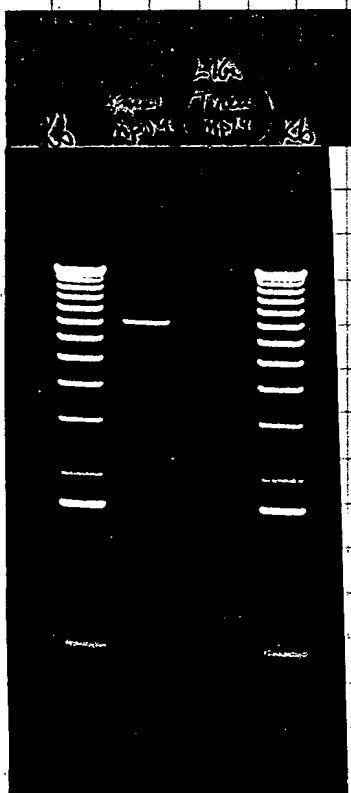

11/21/95

To Page No. _____

Read & Understood by me, <i>May Longo</i>	Date 11/27/95	Invent d by <i>Dr. J. Schmidt</i>	Date 1-26-95
		Recorded by	

56

Project No. 20222


Book No. 3884

TITLE

T. neapolitana 30M

From Page N 3

February 7, 1995 (Tuesday)

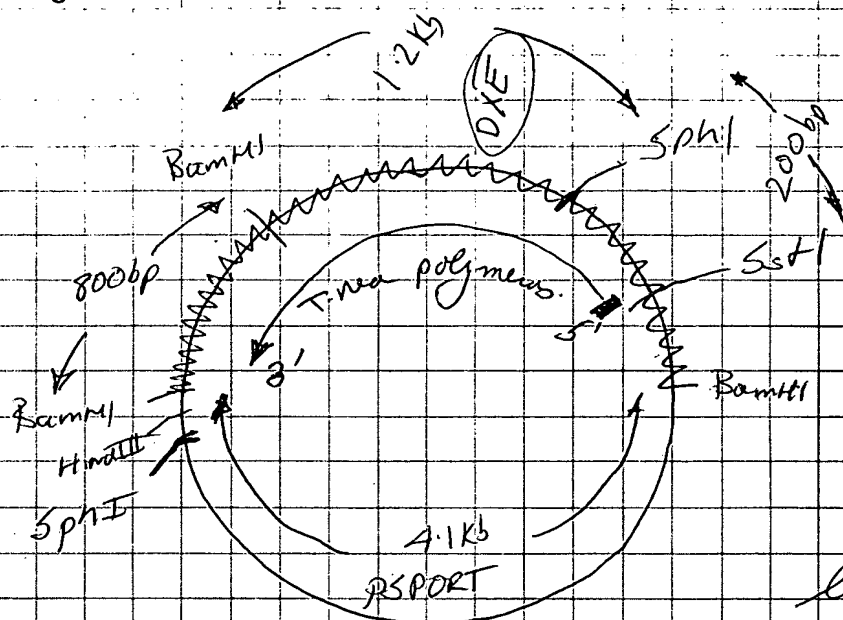

I can clone a SphI fragment into m13mp19 and let m13 determine which direction is best suited.

I can subclone the fragment into an expression vector with SphI / HindIII

DIGEST SCHEME

	T. nea RSPORT	(272 ng/μl) m13mp19	
(Readd) HOH	21 μl ✓	20 μl ✓	9:42 am →
10X Bfr	3 ✓	3 ✓	
DNA	3 ✓	3 ✓	
CMG105 (10 μl) SphI	3 ✓	3 ✓	
(0.1 μl) CAP	0 ✓	1 ✓	9:46 am →
Torax	30 μl	30 μl	

Run on 0.8% Agarose Gel at 75 Volts

To Page

Witnessed & Understood by me,

May Longo

Date

2/16/95

Inv. nted by

Recorded by

[Signature]


Date

2-7-95

g No. 56

February 7, 1995 (Tuesday)

0.8% Agarose Gel (1X TAE); 75 V, 16

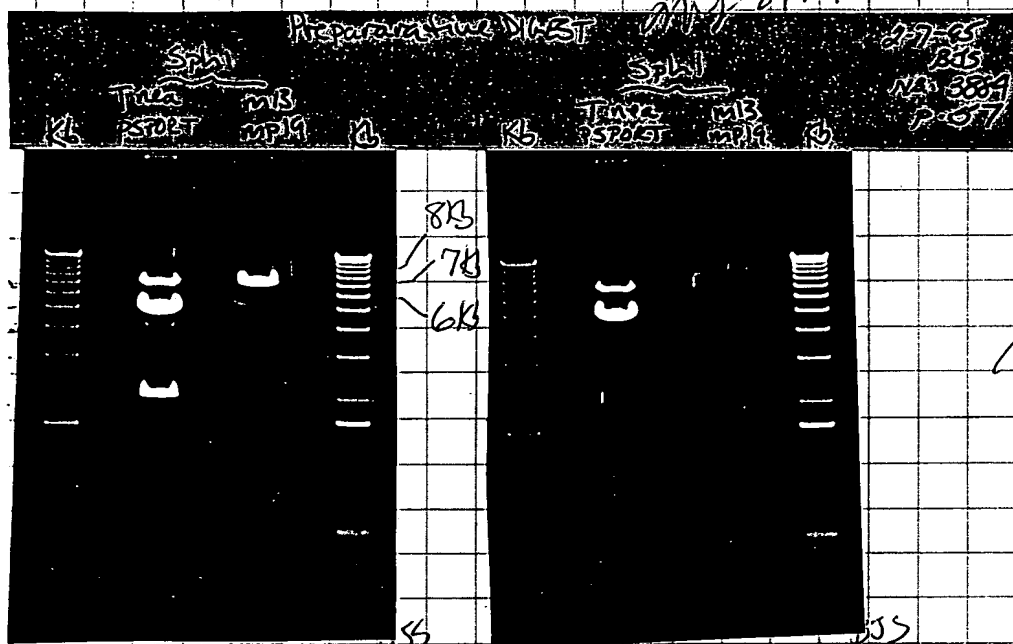

Bands extracted from the gel and the DNA purified away from the agarose using Gene Clean as described by 8.10-10.1. DNA eluted in 14 µl HOH

LIGATION SCHEME

HOH	-	µl	
5X Bfr	4		✓
DNA	14		✓
WAT DNA Ligase	2		✓
TOTAL	20	µl	

Incubated 3:23 pm → 4:03 at room-temperature (~22°C)
3 µl removed for transformation

152 F' IQ Competent Cell Transformation
0.1 µl competent cells + 3 µl ligation (see above)
2 min on ice, 35 seconds at 42°C water bath
1. and 90.1. applied to LB + No Antibiotic plates in 4 ml 0.7% Top Agar + 100 µl 2% X-Gal + 10 µl 100 mM IPTG
incubated 16 hours at 37°C incubator


To Page No. 58

ed & Understood by me, Ney Longo	Date 2/16/95	Invented by [Signature]	Date 2-7-95
-------------------------------------	-----------------	----------------------------	----------------

From Page No. 57

February 8, 1985

- I added 200 μ l of DH5 α F' IQ lawn cells to 10 ml
Archie Brown
→ I added 1 ml of the cells to 8 glass tubes
→ Each tube was inoculated with a clear plaque
and incubated at 37°C (8:00 am →

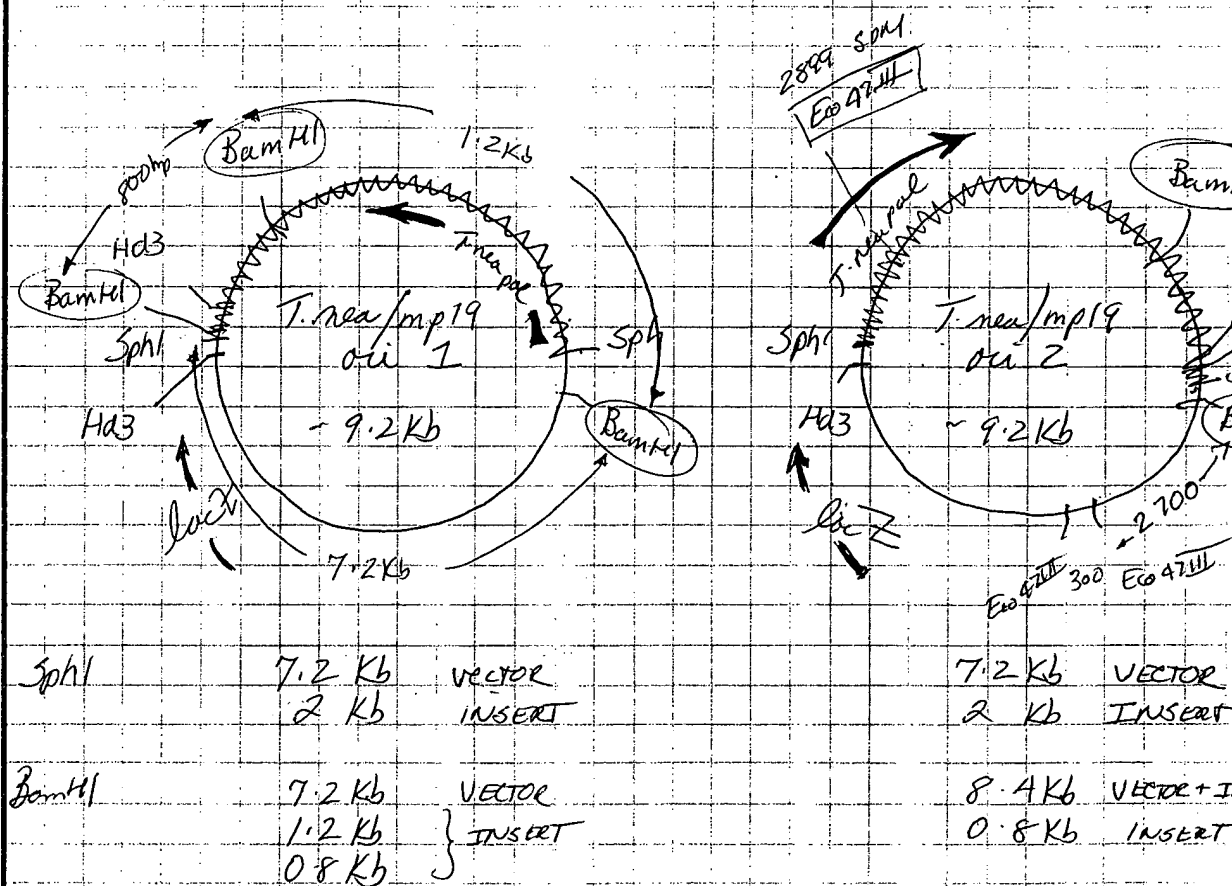

Add *T. nea*/pSPORT as a positive control for both digests

To Page

Witnessed & Understood by me,

M. J. J. J.

Date

2/11/85

Invented by

Recorded by

Date

2-8-85

T-neapolitaner SDM

Book No. 3884

59

ag No. 58

February 8, 1995 (Wednesday)

DIGEST SCHEMES

1.		PER RXN	x 9 =	COCKTAIL		For T-neap/PSPORT
React 6)	HOH	7 μ	x 9 =	63 μ	<input checked="" type="checkbox"/>	control add
	10x Bfr	2	x 9 =	18 μ	<input checked="" type="checkbox"/>	
	DNA	10				T10 E1 20 μ <input checked="" type="checkbox"/>
	(100 μ l) SpH1	1	x 9 =	9 μ	<input checked="" type="checkbox"/>	DNA 2 μ <input checked="" type="checkbox"/>
	Torn	20 μ		90 μ		Torn 22 μ <input checked="" type="checkbox"/>

✓

add 10 μ to reaction

		PER RXN	x 9 =	COCKTAIL	
React 3)	HOH	7 μ	x 9 =	63 μ	<input checked="" type="checkbox"/>
	10x Bfr	2	x 9 =	18	<input checked="" type="checkbox"/>
	DNA	10			
	(100 μ l) BomH1	1	x 9 =	9	<input checked="" type="checkbox"/>
	Torn	20 μ		90 μ	

Continued on page 1 of Notebook 3966

Arani Patel

Gel photo

To Page No. _____

Sed & Understood by me,

Date

Invented by

Date

May Longo

2/16/95

Recorded by


Brenda L. Lamm

2-8-95

Page No. _____

SDM 2899

RF map

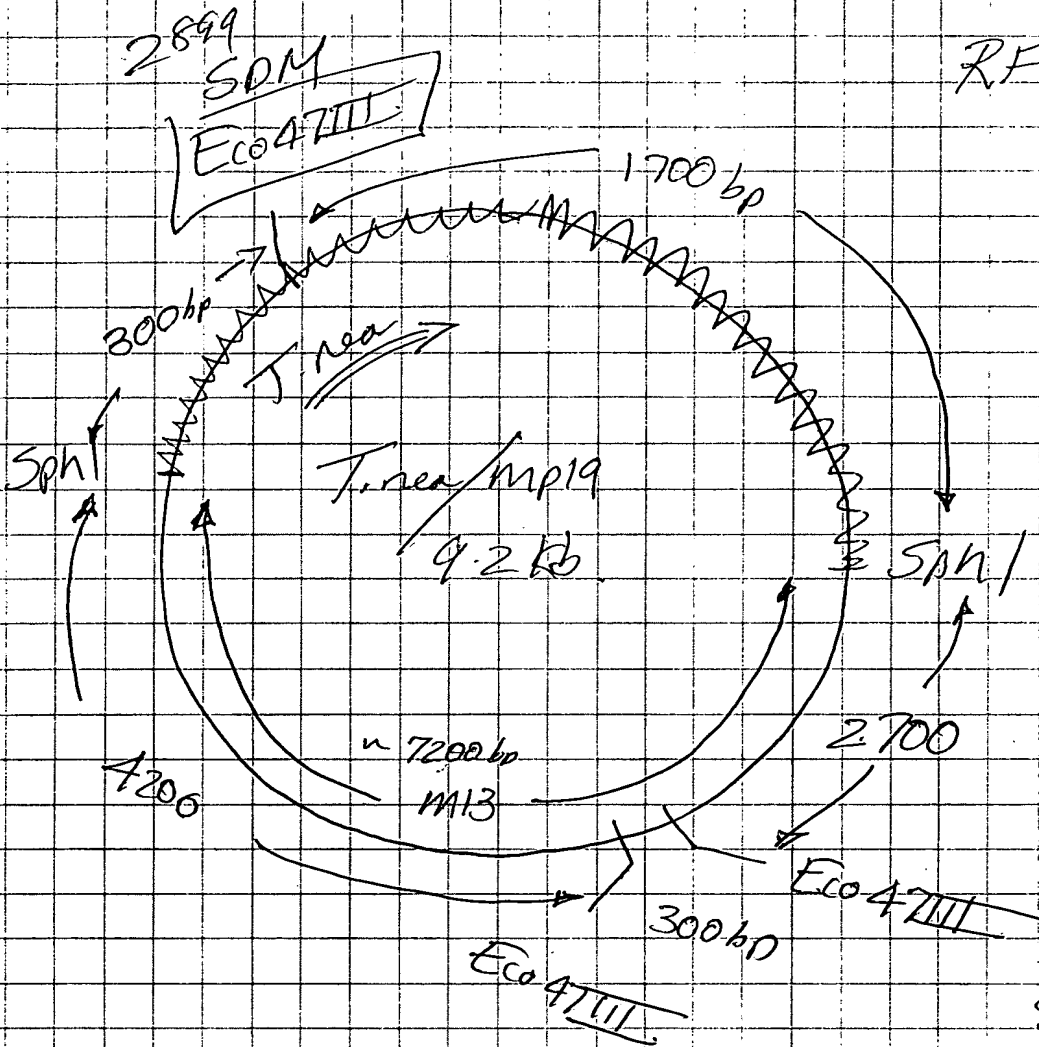

See on Gel

	Eco 47III		1kb	PARENT	MUTANT
PARENT	8.9 Kb	5kb		—	—
	0.3 Kb	4kb			—
		3kb	—		
		2kb	—		
		1.6kb	—		
MUTANT	4.5 Kb	11b	—		
	4.4 Kb	500b			
	0.3 Kb		—		

may be too light to see

Read & Understood by me, Man Jongo	Date 2/16/95	Invented by Dwight. K. Smith	Date 2-16-95
	Page No. _____		

The mutant Phe to Ala

g No. _____

① The same phenyl alanine corresponding to Tag polymerase will be changed to tyrosine

② For exo D will be changed to Alanine (corresponding region of Tag).

Brian cloned the SphI fragment of Tne Pol into M13mp.

I isolated the single stranded DNA from CJ236 as described before in Bio rad manual.

Test 5µl ssDNA

The DNA looks real good.

For D-A (3'-5' exo mutant oligo) is

5' GA | CGT | TTC | AAG | CGC | TAG | GGC | AAA | AGA # 2899
Eco47III site

For Phe → Tyr (O-helix)

~~GA~~ GTA | TAT | TAT | AGA | GTA | GTT | AAC | CAT | CTC | TCC | A
2904

kinased 2899 before.

kinased 2904 as follows:

2µl oligo (210 pmoles)
6µl 5X buffer (350mM Tris pH 7.6, 50mM MgCl₂, 50mM KCl, 5mM P.M.E.)
1µl 10mM ATP
0.5µl T4 Kinase (50)
20.5µl H₂O
5' at 37°C → Heat at 65°C + 3µl TE
T Page 10/22

Sed & Understood by m ,

[Signature]

Date

4/8/95

Invented by

Record d by

[Signature]

Date

3/14/95