

PMC9 / D.V.

Pag N _____

purpose To amplify PMC9 with Deep Vent alone
 tried in Deep Vent buffer only.


10x buffer	100	200 µl	dNTP
dNTP	20	2 µl	Mg
Mg	—	1 µl	primer
primer	10	200 µl	Template
2	10		
Template	4		
420	846	94°	30

30(94°, 30', 56°, 30', 72°, 3'

0	1	2
.001	3	4
.005	5	6
.01	7	8
.05	9	10
.1	11	12
.5	13	14
1	15	16
2	17	18

result: Once again
 Deep Vent alone by itself
 didn't amplify anything
 with this set primers.

See the other set to
 confirm.

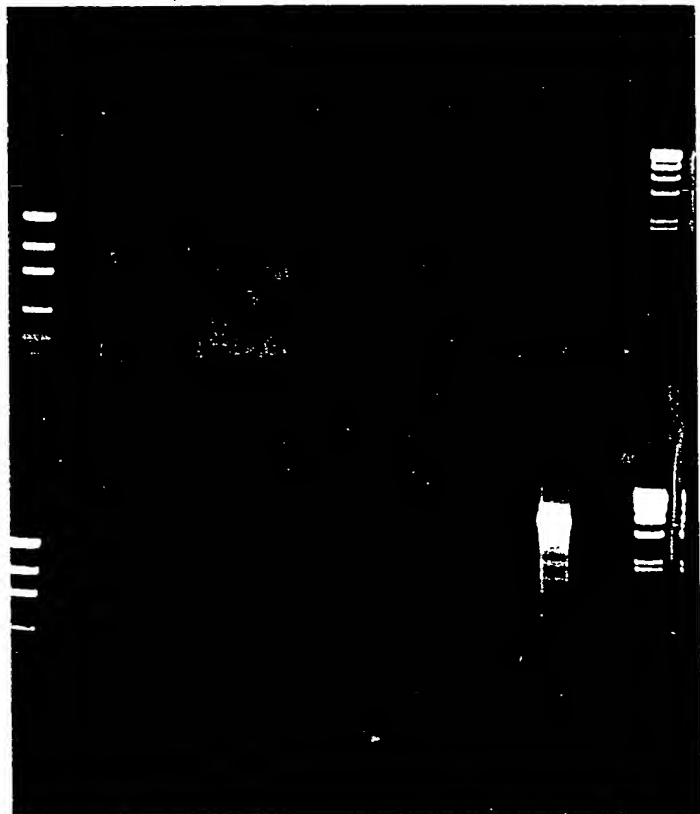

To Page No. _____

Issued & Understood by me,

Date

Invented by

Date

Recorded by

11/22/94

A. Sitarman