


g N. —

Clone 3'→5' exo and F→Y mutation in ptrene as.
 Sst I / Hind III fragment.

Digest ptre 99 with Sst I + Hind III.
 Digest pucTne mutants with Sst I + Hind III + Nde I (to separate almost doublet).

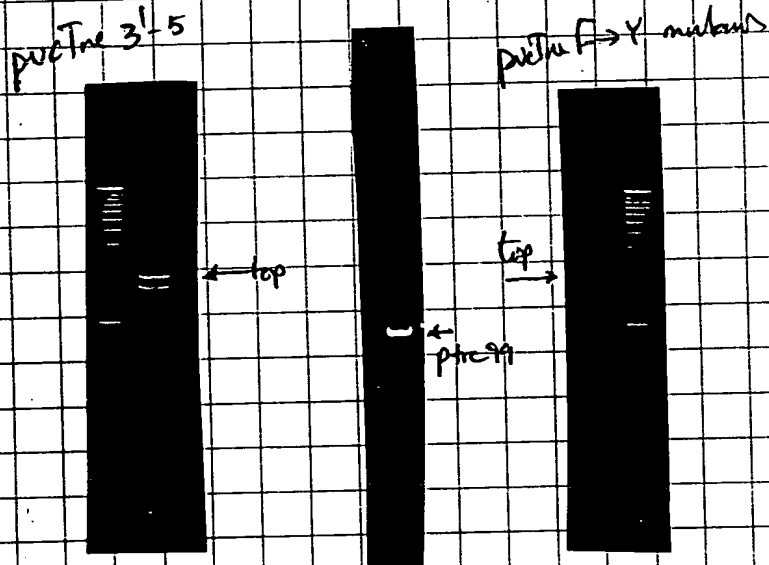

Purify top bands from the mutants and ptre 99.

Ligate → transform D110B.

10% & 90% plating.
 Incubate overnight at 30°C

To Page No. _____

essed & Understood by me,

Date

6/12/95

Invented by

Recorded by

Date

5/16/95