


g N. —

Clone 3'→5' exo and F→Y mutation in pTrcne as.  
 Sst I / Hind III fragment.

Digest pTrc99 with Sst I + Hind III.  
 Digest pucTne mutants with Sst I + Hind III + Nde I (to  
 separate almost doublet).

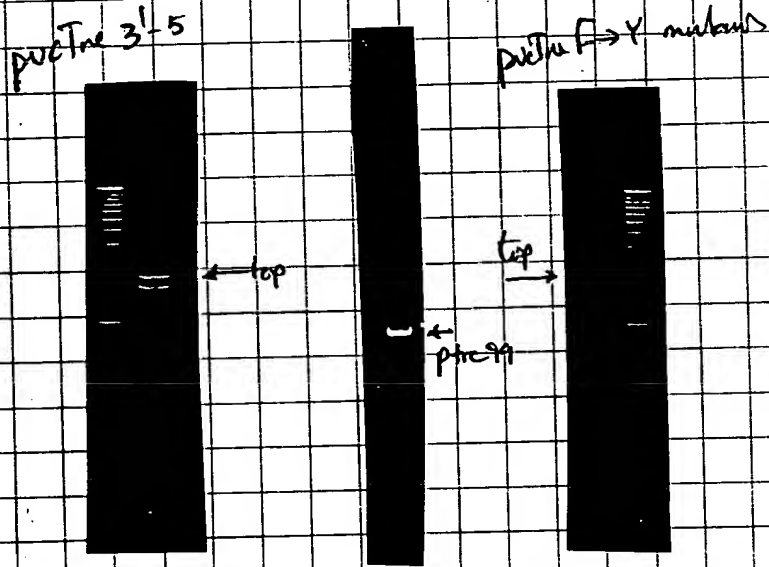

Purify top bands from the  
 mutants and pTrc99.

Ligate → transform D110B.

10% & 90% plating.  
 Incubate overnight at 30°C

To Page No. \_\_\_\_\_

essed & Understood by me,

*[Signature]*

Date

6/12/95

Invented by

*[Signature]*

Date

5/16/95

Recorded by