

The Mutant 3'-5' cko

Project N \_\_\_\_\_  
Book N \_\_\_\_\_

29

Page No. \_\_\_\_\_

Oligo # 2899


D → A mutation.

GA | CGT | TTC | AAG | CGC | TAG | GGC | AAA | AGA

Eco 47III

Kinasing, annealing and Synthesis were done as p. 27.

1 2 3

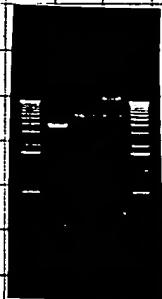

1 = ss DNA (uracil template)

2 = - oligo

3 = + oligo 2899.

Here again a lot of Synthesis without any added oligo.

Nevertheless - proceed with transfection. Dil synthesis reaction 25. Use 2  $\mu$ l for transfection.

10% & 90% plating using same lawn cells.

Inoculate plates at 37°C O/N.

To Page No. \_\_\_\_\_

Read & Understood by me,

Date

Invented by

Date

Recorded by

Deb A. Bhatnagar

4/15/95